


Persee T8DCS Double Beam UV-Vis Spectrophotometer [New]

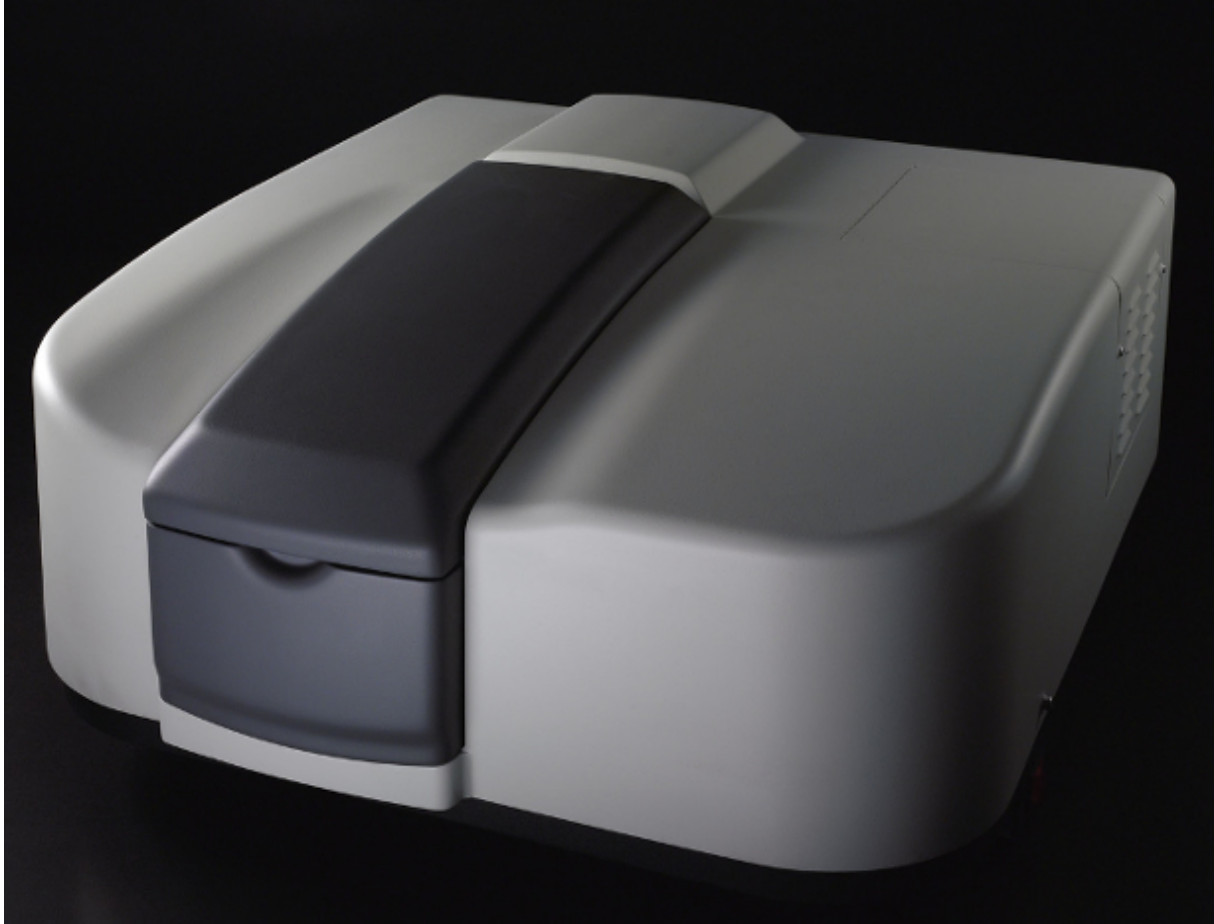

Double Beam
Price with Discount \$7,699.00
Discount

[Ask a question about this product](#)

Manufacturer [Persee](#)

Description

T8DCS is a high performance double beam spectrophotometer with a continuously selectable spectral bandwidth from 0.1-5nm.

The Czerny-Turner Monochromator with a holographic grating keeps stray light to a minimum and

offers excellent optical resolution.

- The use of a photomultiplier tube as detector delivers exceptional sensitivity.
 - The T8DCS' true double-beam optical system coupled with an efficient and well proven electronic control system ensures high stability and low background noise.
 - User friendly design allows easy light source replacement and routine maintenance.
 - Large sample compartment design enables use of a wide range of accessories.
 - Wavelength: 190-900nm
 - Variable Bandwidth: 0.1, 0.2, 0.5, 1, 2, 5nm
-